


INSTRUCTIONS FOR USE
REDEMPLO [ree-DEM-plo]
(plozasiran)
injection, for subcutaneous use
single-dose prefilled syringe
25 mg/0.5 mL

This Instructions for Use contains information on how to inject REDEMPLO (plozasiran).

Read this Instructions for Use before you start using your REDEMPLO prefilled syringe and each time you get a refill. There may be new information. This information does not take the place of talking to your healthcare provider about your medical condition or your treatment. Your healthcare provider should show you or your caregiver how to use your REDEMPLO prefilled syringe the right way. Call your healthcare provider or pharmacist if you have any questions.

Guide to Parts

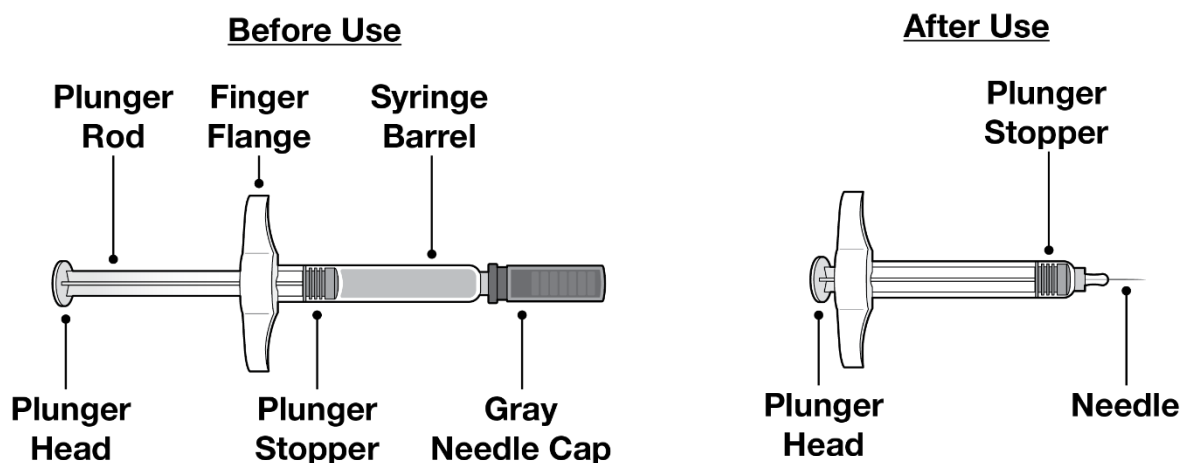

Each carton contains a single-dose REDEMPLO (plozasiran) Injection prefilled syringe for single use under the skin (subcutaneous) injection.

Important Information You Need to Know Before Injecting REDEMPLO

- **REDEMPLO is for subcutaneous injection only**
- Each carton contains a 1-time (single-dose) REDEMPLO (plozasiran) prefilled syringe.
- **Do not** use your REDEMPLO single-dose prefilled syringe if the carton appears damaged or the tamper evident seal on the carton is broken.
- **Do not** remove the gray needle cap until you are ready to inject (Step 6).
- **Do not** use if the single-dose prefilled syringe appears damaged.
- Throw away (dispose of) your REDEMPLO prefilled syringe in a sharps disposal container right away after use (Step 10: Disposing of REDEMPLO).

Storing REDEMPLO

- Keep the REDEMPLO prefilled syringe in the original carton until ready to use.
- Store the REDEMPLO prefilled syringe in the refrigerator between 36°F to 46°F (2°C to 8°C) in the original carton.
- REDEMPLO prefilled syringe can also be stored at room temperature between 68°F to 77°F (20°C to 25°C) in the original carton for up to **30 days**.
- Do not let REDEMPLO reach temperatures above 77°F (25°C).
- **Throw away** the REDEMPLO prefilled syringe if kept at room temperature and not used within 30 days.

If the product is not stored in any of the above conditions, throw away the prefilled syringe in a sharps disposal container and use a new prefilled syringe.

Keep REDEMPLO prefilled syringe and all medicine out of reach of children.

Preparing to Inject REDEMPLO

Step 1: Gather all materials needed for your REDEMPLO injection




**Prefilled syringe
inside carton**



**Alcohol
wipes**



**Cotton ball
or gauze pad**

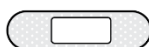

**Adhesive
bandage**



**Sharps disposal
container**

Figure A

On a clean, well-lit, flat work surface, place:

Provided in the carton (see **Figure A**):

- 1 REDEMPLO prefilled syringe

Other materials not provided (see **Figure A**):

- Alcohol wipes
- Cotton ball or gauze pad
- Adhesive bandage
- Sharps disposal container

Step 2: Prepare to use REDEMPLO prefilled syringe

Grasp by syringe
barrel to remove

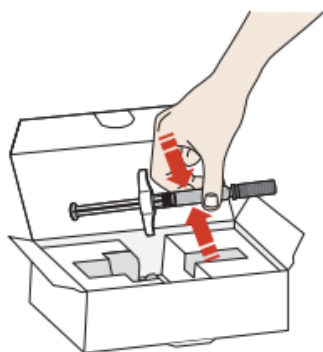


Figure B



Expiration Date

Figure C

Wait to reach
room temperature



Figure D

- If refrigerated, remove the carton from the refrigerator.
- Open the carton lid and remove the syringe by the syringe barrel and place on a flat surface (see **Figure B**).
 - **Do not** use the prefilled syringe if tamper evident seal on the carton is broken.
 - **Do not** pick up or pull the prefilled syringe by the plunger rod or gray needle cap.
- **Check the expiration date** on the REDEMPLO prefilled syringe (see **Figure C**).
 - **Do not** use if the expiration date has passed.
- Wait 30 minutes for the prefilled syringe to reach room temperature before injecting (see **Figure D**).
 - **Do not** try to warm the prefilled syringe by using a heat source such as hot water or a microwave.
 - **Do not** remove the gray needle cap from the prefilled syringe until you are ready to inject.

Step 3: Check the medicine and syringe



Figure E

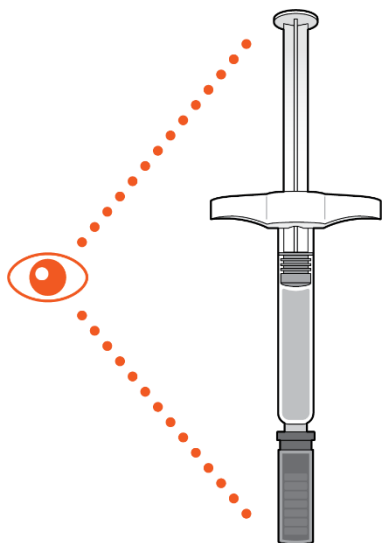


Figure F

- Check the REDEMPLO medicine in the prefilled syringe (see **Figure E**).
 - The medicine should be clear and colorless to yellow.
 - **Do not** use the prefilled syringe if the medicine is cloudy, discolored or contains particles.
 - It is normal to see air bubbles in the medicine.
- Inspect the prefilled syringe (see **Figure F**).
 - **Do not** use the prefilled syringe if any part appears cracked or broken.
 - **Do not** use the prefilled syringe if the gray needle cap is missing or not securely attached.
 - **Do not** use the prefilled syringe if it has been dropped onto a hard surface because the syringe may be damaged.
- In any of the above cases, call 1-844-REDEMPLO or visit www.redemplo.com.

Injecting REDEMPLO

Step 4: Choose your injection site

Upper-outer arm for healthcare providers and caregiver only

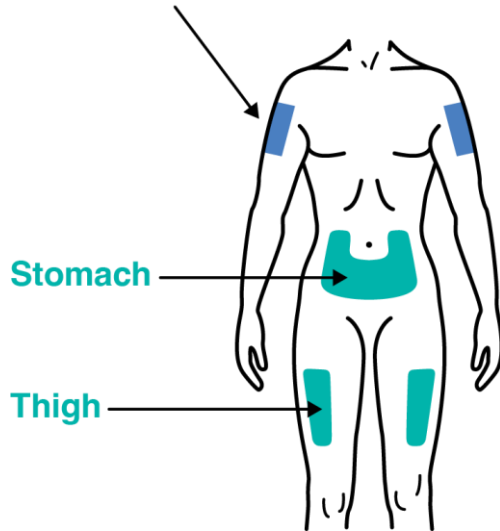


Figure G

- **Choose an injection site on the (see Figure G):**

- Upper leg (thigh) or
- Stomach (abdomen) except for a 2-inch area right around the belly button (navel).

Only your healthcare providers or caregivers can use the outer area of the upper arm (see Figure G).

- **Do not inject:**

- within 2 inches of the belly button (navel).
- into skin that is damaged (tender, bruised, red, hard, or cut).
- into areas with scars, stretch marks.

Step 5: Wash hands and clean the injection site

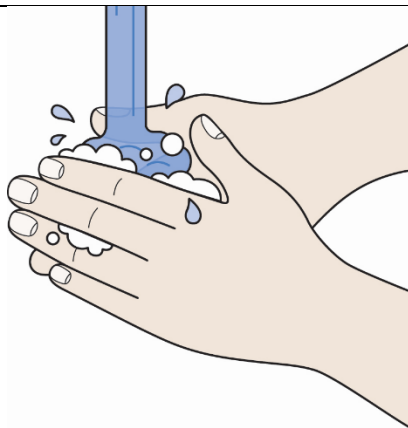


Figure H

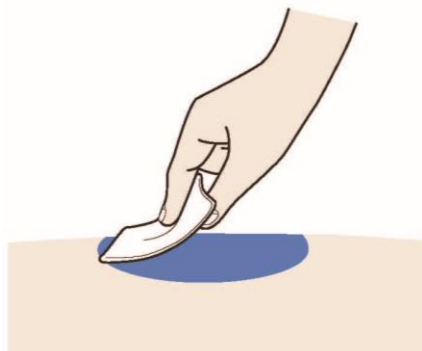


Figure I

- Wash your hands thoroughly with soap and water (see **Figure H**).

- Clean your injection site with an alcohol wipe. Let your skin dry before injecting (see **Figure I**).
 - **Do not** touch the cleaned skin before injecting.

Step 6: Remove the needle cap

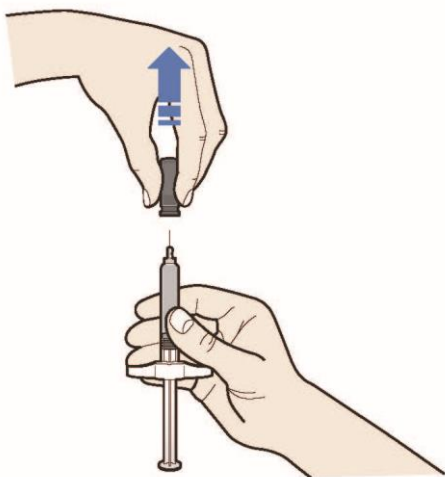


Figure J

- Hold the prefilled syringe by the syringe barrel, with the needle facing away from you.
- Pull the gray needle cap straight up and away from your body. (see **Figure J**).
 - **Do not** twist or bend the gray needle cap.
 - Avoid pushing the plunger head before you are ready to inject.
 - **Do not** let the needle touch any surface.
 - **Do not** put the needle cap back onto your prefilled syringe.

Step 7: Pinch the skin and Insert the needle

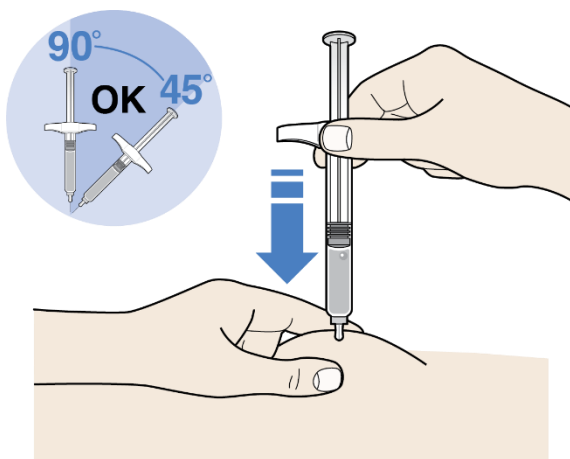


Figure K

- Hold your prefilled syringe by the finger flange in one hand.
- Gently pinch and hold a fold of skin at the injection site with your other hand.
- Insert the needle at a 45°(degree) to 90° angle (see **Figure K**).
- **Keep the skin pinched while you are inserting the needle and during injection.**
 - **Do not** place your finger on the plunger head before the injection.

Step 8: Hold skin pinch and Push down the plunger head

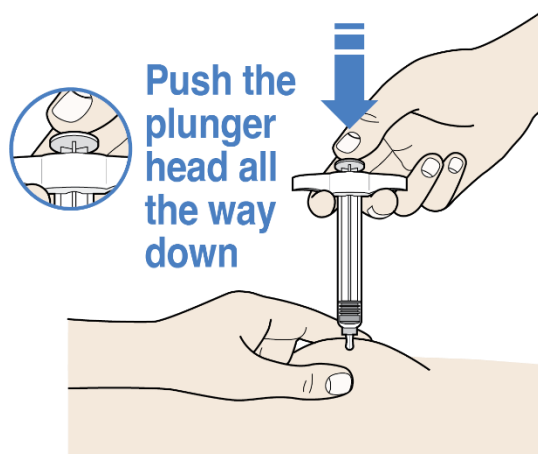


Figure L

- While pinching the skin, push the plunger head all the way down using slow and constant pressure (see **Figure L**).

Step 9: Complete the Injection

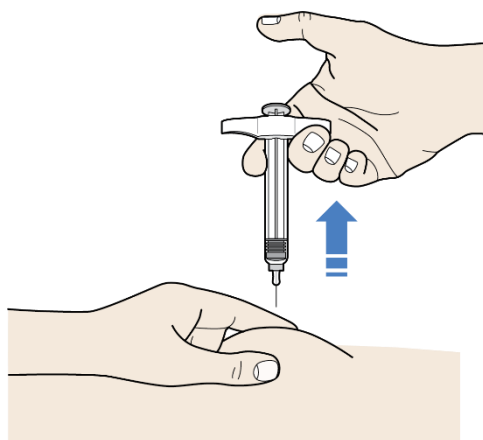


Figure M

- After the plunger head is pushed all the way down, take the needle out of the skin, by gently lifting the syringe up and off the skin until the needle is completely removed (see **Figure M**).
 - **Do not** pull the plunger head up. Lift the whole syringe straight up.
 - **Do not** rub the injection site.
- There may be a small amount of blood or liquid where you injected. This is normal.
 - If needed, press a cotton ball or gauze pad on the area and apply an adhesive bandage.
- Place the used cap and syringe in a sharps disposal container right away (see section – “Disposing of REDEMPLO”).
 - **Do not** put the needle cap back onto the syringe.

Step 10: Disposing of REDEMPLO



Figure N

Do not use any medicine that is left in the used syringe.

- Put your used syringes and needle cap in an FDA-cleared sharps disposal container right away after use.
- **Do not** throw away (dispose of) loose needles and syringes in your household trash (see **Figure N**).
- If you do not have an FDA-cleared sharps disposal container, you may use a household container that is:
 - made of a heavy-duty plastic,
 - can be closed with a tight-fitting, puncture-resistant lid, without sharps being able to come out,
 - upright and stable during use,
 - leak-resistant, and
 - properly labeled to warn of hazardous waste inside the container.
- When your sharps disposal container is almost full, you will need to follow your community guidelines for the right way to dispose of your sharps disposal container. There may be state or local laws about how you should throw away used needles and syringes.
- For more information about safe sharps disposal, and for specific information about sharps disposal in the state that you live in, go to the FDA's website at:
<http://www.fda.gov/safesharpsdisposal>.

Do not dispose of your used sharps disposal container in your household trash unless your community guidelines permit this. **Do not** recycle your used sharps disposal container.

For more information, call 1-844-REDEMPLO or visit www.redemplo.com.

Distributed by:
Arrowhead Pharmaceuticals, Inc.
Pasadena, CA 91105
Phone Number: 1-844-REDEMPLO

This Instructions for Use has been approved by the U.S. Food and Drug Administration.

Approved: 11/2025